


Digital Crockmeter ATL:0019

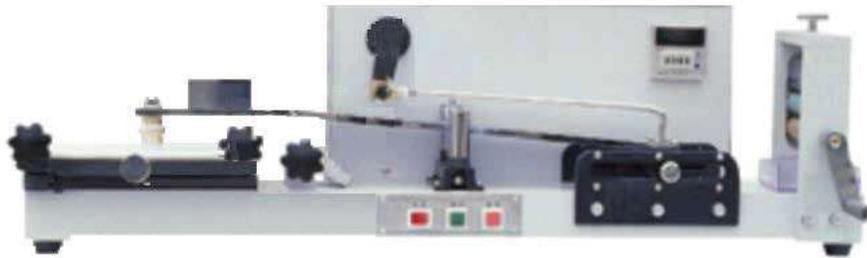

The specimen is clamped to the baseplate, where a peg is moved back and forth over the test specimen to determine colorfastness from dry or wet tests equipped with hand crank.

Applications:

- Textiles
- Carpets

Features:

Testing Times:	- 1-9999 (LED)
Rubbing Travel:	- 104mm ± 3mm
Rubbing Pressure:	- 9N
Rubbing Speed:	- 60r/min
Rubbing Head:	- Round 16mm Ø Rectangle 19mm x 25mm

Standards:

- AS 2001.4.3
- AS / NZ 2111.19.1
- AATCC - 8-(1972), AATCC 165
- BS 1006
- GB/T3920
- ISO 105-D02, ISO105-X12,
M&S C8, M&S C8A, M&S C25, M&S C52
- IWS TM165, JIS K6328

Next 10

Benefits:

- Easy to use
- Fast results
- Accurate

Connections:

- AC 220V

Dimensions:

- H: 280mm
- L: 770mm
- D: 250mm
- Weight: 25kg

Accessories:

- Steel specimen holder
- 17mm test finger.

Regulatory compliances



Corporate Social Responsibility

Analytical Foundation is a Nonprofit Organization (NGO) found for the purpose of:



**Analytical
Foundation**

1. Research & Innovation Scientist's awards / QC Professional Award : Quality life is possible by innovation only and the innovation is possible by research only, hence ANALYTICAL FOUNDATION is committed to identify such personalities for their contributions across various field of Science and Technology and awarding them yearly. To participate for award, send us your details of research / testing / publication at info@analyticalfoundation.org

2. Improving quality of life by offering YOGA Training courses, Work shops / Seminars etc.

3. ANALYTICAL FOUNDATION aims to DETOXYFY human minds, souls and body by means of Yoga, Meditation, Ayurveda, Health Care, Awards, Media, Events, Camps etc.

Reach us @



Corporate & Regd. Office: Analytical House, # E67 & E68, Ravi Park, Vasna Road, Baroda, Gujarat 390 015. INDIA	HPLC Solutions T: +91 265 2253620 +91 265 2252839 +91 265 2252370 F: +91 265 2254395	MultipleLabs	Analytical Bio-Med E: info@hplctechnologies.com info@multiplelabs.com info@analyticalgroup.net info@analyticalbiomed.com	Analytical Distributors W: www.ais-india.com www.analycalgroup.net www.hplctechnologies.com www.multiplelabs.com	Analytical Foundation (Trust) Sales & Support Offices: across the country : Distributors & Channel partners World Wide
--	---	---------------------	--	---	---